




OXYGEN UNLIMITED

NEO 5

New Generation, High-efficiency
Cost-effective, Reliable

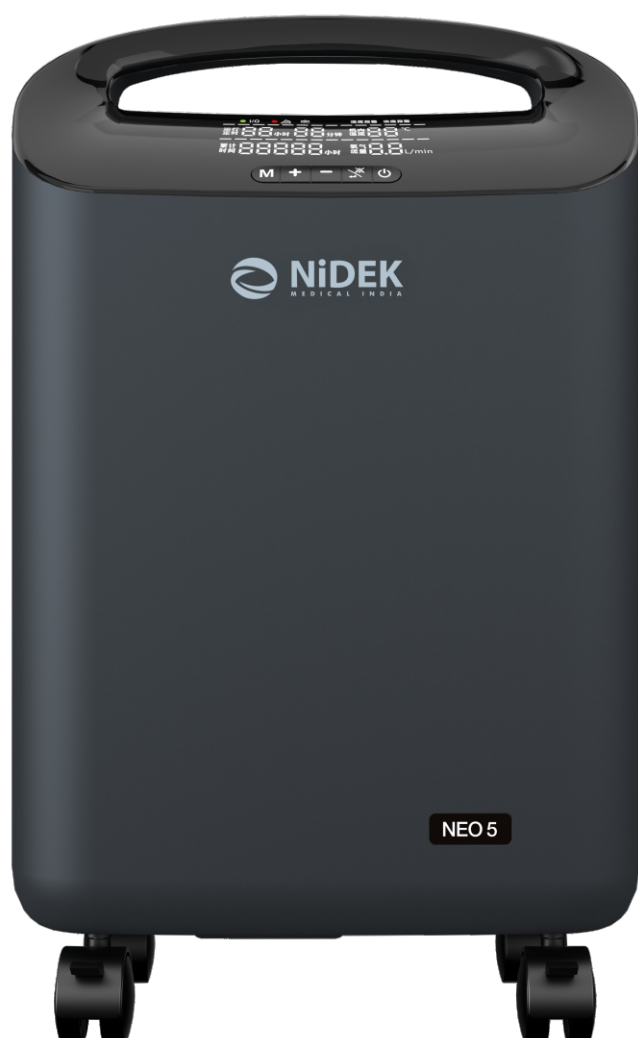

The NEO 5 is the newest and most advanced evolution in Nidek India's Oxygen Concentrator line. Designed with innovation at its core, it's now lighter, more compact, and more capable than ever before.

Superior quality features:

- **Lite:** Only 14.4 Kg
- **Quiet:** Ultra quiet operation, sound level ≤ 50 dBA
- **Compact:** $32.2 \times 22.2 \times 52.2$ cm
- **Accurate:** Electronically controlled flow rate
- **Economical:** Very low power consumption only 300 W

Other features:

- Rugged, reinforced construction with durable castors
- High-efficiency compressor
- Easy to replace air inlet filter and cabinet filter
- It can be easily adjusted to the prescribed flow rate with a single touch of button
- Full purity in less than 3-minutes
- Offers significant reduction in preventive and curative maintenance time
- Easy usage and maintenance

Technical Parameters

Capacity	0 to 5 LPM
Continuous Flow Settings	1.0 to 5 LPM
Size (in mm)	322 W × 222 D × 522 H
Oxygen Concentration	93 ± 3%
Oxygen Output Pressure	55 ± 10% Kpa
Sound Level	≤ 50 dBA
Power Specifications	230 VAC, 50 Hz
Power Consumption	300 Watts
Filters	Cabinet, Compressor Intake and Bacteria
Weight	14.4 Kg
Nebulizer Function	Yes
Oxygen Purity Indicator	Yes
Alarm	Low Purity Alarm Over Heat Alarm Power Failure Alarm



Trusted Name Of Nidek



Accurate Oxygen Monitoring System



High Quality Sieve For Superior Performance



Nationwide Service Network



Nidek Medical India Private Limited

Head Office: 238 B, A.J.C. Bose Road, Kolkata 700 020

Tel: (033) 4017 7000, E-mail: sales@nidekindia.com

Website: www.nidekindia.com

Branch Locations:

New Delhi | Mumbai | Chennai | Hyderabad | Lucknow

Pune | Nagpur | Ahmedabad | Chandigarh | Indore

Bangalore | Kochi | Guwahati

Customer Care: +91 93395 66444